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# CONTENTS

Safety precautions	31
Introduction	32
Precautions before use	33
Functions	33
Usage	34
Practical considerations when taking a temperature	34
How to take a temperature	35
Important information when taking a temperature	35
Description	36
Setting and function	37
Instructions for calibration	39
Changing the batteries	39
Technical characteristics	40
Troubleshooting	41
Warranty	43

## SAFETY PRECAUTIONS

- Follow the maintenance advice stipulated in this instruction manual.
- This device may be used for professional purposes or for personal home use.
- This device must only be used in an ambient temperature range of between 10°C and 40°C.
- This device must always be kept in a clean, dry area.
- Do not expose this thermometer to electric shocks.
- Do not expose this thermometer to extreme temperature conditions of  $> 55^{\circ}\text{C}$  or  $< -20^{\circ}\text{C}$ .
- Do not use this device in relative humidity higher than 85%.
- The protective glass over the lens is the most fragile part of the thermometer.
- Do not touch the glass of the infrared lens with your fingers.
- Clean the glass with a cotton bud lightly moistened with 95% alcohol.
- Do not expose the thermometer to sunlight or to water.
- Never drop the device.
- Do not attempt to repair this device yourself. If you have any problem, please, contact authorised service department.

The manufacturer reserves the right to alter the specification of the product without prior notification.

## INTRODUCTION

The non-contact infrared thermometer Helpmation has been developed by using the latest infrared technology. This technology allows temporal artery (TA) temperature to be taken at a distance of max. 3 cm away from the forehead. It has been demonstrated that this method of TA temperature measurement is more precise than the tympanic thermometry and better tolerated than rectal thermometry (1). However, as with other types of thermometer, it is essential to use the thermometer properly in order to obtain reliable and stable results.

(1) „Greenes D, Fleischer G. Accuracy of a Noninvasive Temporal Artery Thermometer for Use in Infants. Arch Pediatr Adolesc Med 2001, 155:376.“

## NORMAL TEMPERATURES ACCORDING TO AGE

Age	°C	°F
0-2 years	36.4-38.0	97.5-100.4
3-10 years	36.1-37.8	97.0-100.0
11-65 years	35.9-37.6	96.6-99.7
> 65 years	35.8-37.5	96.4-99.5

## PRECAUTIONS BEFORE USE

The thermometer Helpmation is pre-set at the factory. It is not necessary to calibrate the device when starting.

In order to obtain reliable and stable results, you are advised each time there is a significant change in the ambient temperature due to a change in environment, to allow the thermometer to acclimatize to this ambient temperature for 15 to 20 minutes before using it. It is important to allow 1 minute interval between two measurements.

## FUNCTIONS

1. Special design to take the Human Body Temperature from a distance.
2. Reliable and stable measurement, thanks to the advantage Infrared Detection System.
3. Audible alarm if temperature is too high.
4. Memorize the last 32 temperature measurements.
5. Backlit LCD display.
6. Temperature can be displayed in either Celsius or Fahrenheit.
7. Automatic power-off (30 seconds) to save energy.
8. Practical, easy to use.
9. Longevity use (40.000 reading).

The thermometer can also be used to measure the temperature of a baby-bottle or bath, or room temperature (by using the Surface Temp or Room function). This function is in accordance with the Directive 89/339/EEC Electromagnetic compatibility.

## USAGE

1. Install the battery (please see Changing the battery)
2. For the first use or when inserting new battery wait around 10 minutes for the warm-up of the device.
3. Aim towards the forehead, from a distance of max. 3 cm, press the measuring button SCAN, the temperature is displayed in 1 second.
4. Before taking the temperature, make sure to remove hair and perspiration from the forehead.

## PRACTICAL CONSIDERATIONS WHEN TAKING A TEMPERATURE


In order to ensure that precise and accurate temperature measurements are obtained, it is essential that each user has received adequate information. Temperature should be taken in a neutral context. The patient must not have undertaken vigorous physical activity prior to taking the temperature.

Be aware of physiological variations in temperature which must be taken into consideration when evaluating the results:

- temperature increases by  $0,5^{\circ}\text{C}$  between 6am and 3pm
- women have a temperature that is higher, on average, by around  $0,2^{\circ}\text{C}$
- women temperature also varies in accordance with their ovarian cycle. It rises by  $0,5^{\circ}\text{C}$  in the second half of the cycle and at the early stages of pregnancy.
- when sitting, temperature is lower by about  $0,3 - 0,4^{\circ}\text{C}$  than when standing

## HOW TO TAKE A TEMPERATURE

Aim at the forehead, over the right temporal region, from a distance of max. 3 cm, press the button SCAN and the temperature is instantly displayed.

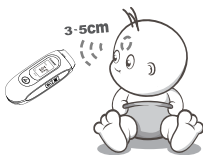The reliability of the measurement cannot be guaranteed if the temperature is measured over another part of the body (arm, torso).



## IMPORTANT INFORMATION WHEN TAKING A TEMPERATURE

NE

Please observe the following before any temperature measurement to ensure a stable and reliable result:

- push back hair from the forehead
- wipe away any perspiration and make-up from the forehead
- avoid any drafts (e.g. from nasal specs, air conditioning ...)
- allow a 1 minute interval between two measurements
- each time there is a significant change in the ambient temperature due to a change in environment, to allow the thermometer to acclimatise to this ambient temperature for at least 15 minutes before using it

# DESCRIPTION
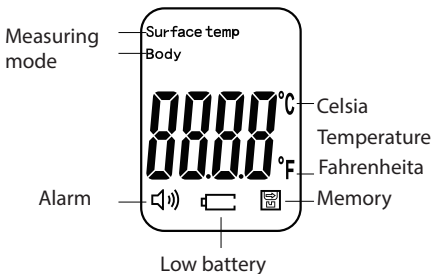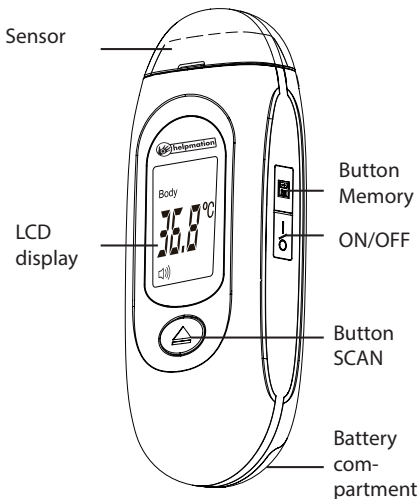
# SETTING AND FUNCTION

## SWITCH ON THE THERMOMETER

If the thermometer is off, press the button ON/OFF and the thermometer switches on (LCD display lights up). Press SCAN measure the desired temperature. Default setting of the thermometer is „BODY“ mode and measuring in °C.

## Setting and function of menu

- **Choosing the temperature °C or °F** - in the switch off state, press the ON/OFF button for two seconds and the screen will display °C or °F. Select „Memory“ button for °C and select „Scan“ button for °F. Confirm by press ON/OFF button. And press ON/OFF button three times to exit the settings.
- **Buzzer ON/OFF** - in the switch off state, press the ON/OFF button for two seconds and the screen will display °C or °F. Press ON/OFF button one time to get Buzzer ON/OFF. Select „Scan“ button the sound will open and display „ON“ and „Alarm“ ikon, select „Memory“ button the sound will close and display „OFF“. Confirm by press ON/OFF button. And press ON/OFF button two times to exit the settings.

- **Mode choice BODY/SURFACE** - in the switch off state, press the ON/OFF button for two seconds and the screen will display °C or °F. Press ON/OFF button two times to get mode choice. Select „Scan“ button to choose SURFACE MODE and select „Memory“ button to choose BODY MODE. Confirma by press ON/OFF button. And press ON/OFF button one time to exit the settings.

## IMPORTANT!

**To obtain the internal temperature of the body, always use the „BODY“ mode. Please, make sure to select „Sufrage temp“ mode for an external area reading.**

### **Data memory**

In the switch on state, press „Memory“ button, which will display the last temperature, and allows for a view of the last 32 measurements taken.

### **Delete the memory**

In the switch on state, press „Memory“ button for 4 seconds, will delete all memory data. Then press „Memory“ button again will display „CLr“.


## CALIBRATION OF THE THERMOMETER

Thermometer factory-set. Depending on various skin types and thickness, there may be temperature differences..


In the switchoff state, press ON/OFF button for two seconds and the screen will display °C or °F. Press ON/OFF button three times to get recalibration set. Select „Scan“ button to add per 0.1 °C or °F. Confirma by press ON/OFF button and exit the settings.

## CHANGING THE BATTERY

When the LCD screen displays „Battery icon“, the battery should be replaced for the new one. Open the lid and change the battery, taking great care with the correct positioning. A mistake with this could cause damage to the unit and compromise the guarantee of your thermometer.

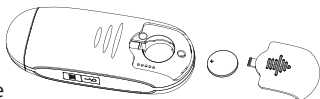Use the battery cover (or toothpick) to press the metal shrapnel in the battery holder, the battery will pop up automatically.
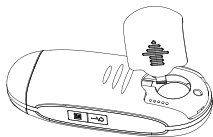
# TECHNICAL CHARACTERISTICS

Normal using condition:

Ambient temperature: 10°C - 40°C (50°F - 104 °F)

Relative humidity: < 85%

Model	RC009
Dimensions	125 x 49 x 27 mm
Weight	64 g (without battery)
Power	DC3 V (1 pc CR2032)
Consumption	< 300 mW
Measuring distance	max. 3 cm
Automatic power-off	30 sec.

Measuring range:

In Body mode 32°C - 42,9 °C (89,6 °F - 109,2 °F)

In Surface Temp mode 0 °C - 60 °C (32 °F - 140 °F)

In Room mode 0 °C - 60 °C (32 °F - 140 °F)

Precision:

36°C ~ 39°C (96,8 °F ~ 102,2 °F) =  $\pm 0.2^{\circ}\text{C}$  ( $\pm 0,4^{\circ}\text{F}$ )

According to ASTM Standard E1965-1998 (2009)

## Precision

From 34°C to 35,9°C = $\pm 0,3^{\circ}\text{C}$ From 93,2°F to 96,6 °F = $\pm 0,6^{\circ}\text{F}$	According to ASTM Standard E1965-1998 (2009)
From 36°C to 39°C = $\pm 0,2^{\circ}\text{C}$ From 96,8°F to 102,2 °F = $\pm 0,4^{\circ}\text{F}$	
From 39°C to 42,5°C = $\pm 0,3^{\circ}\text{C}$ From 102,2°F to 108,5 °F = $\pm 0,6^{\circ}\text{F}$	

The thermometer takes temperature readings below 32 °C or above 42,9 °C (89,6 °F to 109,2 °F) but precision is not guaranteed outside of this range.

## TROUBLESHOOTING

If you have problems while using your thermometer, please refer to this guide to help resolve the problem. If the problem persists, please contact our authorised service department.

- **The screen displays temperature higher than 42,9 °C (109,2°F):**  
The temperature is in Fahrenheit. Change the measurement to Celsius by pressing °C/°F button.
- **The screen displays temperature lower than 32°C (89,6°F):**  
To take the surface temperature, press the "BODY" button and set to the reading called "Body", If the device is in Surface Temp Mode, the 32°C (89.6°F) temperature displayed is showing the external temperature of your body, rather than the internal.
- **The screen displays the message „HI“:**  
When using the Thermometer, the message "HI" can show on the screen. In this case, the temperature is above the measurement range selected, either above 42.9 °C (109.2°F) in Body Mode.

- **The screen displays the message „Lo“:**

When using the RC006 Thermometer, the message “LO” can show on the screen. In this case, the temperature analyzed is under the measuring range selected, either less than 32°C (89.6°F) in Body Mode.

This message displays for various reasons:

- push back hair and wipe the make-up from the forehead
- measurement distance is longer than 3 cm
- allow a 1 min. interval between two measurements
- avoid any drafts (nasal specs, air conditioning)
- make sure that you measure the desired temp. in the proper setting of the thermometer (BODY, SURFACE)

### Contents:

- 1 pc non-contact thermometer
- 1 pc CR2032 battery
- 1 pc instruction manual

This product is in compliance with the essential requirements and other relevant provisions with EU Directives 93/42/EEC.

Do not dispose the product in household waste at the end of its useful life. Please, dispose of this device at your applicable collection point for the recycling of electrical & electronic equipments waste.

Guangzhou Jinxinbao Electronic Co., Ltd.  
Address: No. 38 Huanzhen Xi Road, Dagang Town, Nansha, Guangzhou, China

Wellkang Ltd.,  
Address: Suite B, 29 Harley Street London, W1G 9QR, England, United Kingdom

# WARRANTY

Model:

Serial number:

Date of purchase:

Signature and stamp dealer:

NE

The product is guaranteed for 24 months.

The warranty covers defects arising as a result demonstrably defective material, faulty design or bad design.

The warranty does not cover defects caused by transport, normal use of the product (in particular the loss of battery capacity), failure to follow the instructions indicated in the instructions, poor maintenance, mechanical damage, unauthorized intervention in the appliance service. Repairs may only be performed by an authorized service center.

Specifications are based on a condition known at the time of publication.

Technical alterations.

BLAKAR trading s.r.o.  
Czech Republic  
[www.helpmation.com](http://www.helpmation.com)